


BLOOD PIPETTE SHAKER

PS-600T

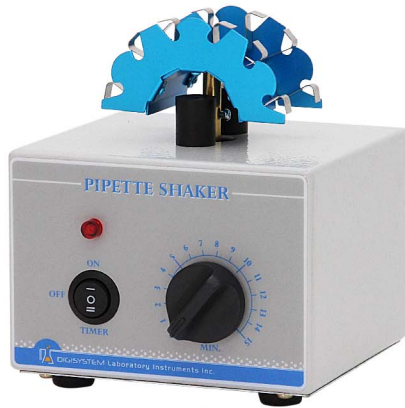

PS-600T

Features:

- Blood pipette shaker is an electric vibrator for blood diluting pipettes which provides an uniform mixing of the chemicals with blood to be tested due to its controlled agitation.
- PS-600T has a pipette holder made of stainless steel spring clips.
- PS-600T can provide continuous or timing operation (0-15min.)
- Quality certificate: ISO 9001, ISO 13485.

Specifications:

Model	Capacity	Switch	Fixed speed	Weight (kg) (N.W./G.W.)	Overall dimension (mm) (WxDxH)
PS-600T	6 pipettes	3 positions (ON/OFF/0-15 min. timer)	2,500rpm	1.8 / 2.2	150x135x160

Power supply: AC 110V 60Hz (220V 50/60Hz available)

SUSPENSION ROTATOR

SR-100



SR-100

Features:

- 3D gyratory movement provides consistently uniform specimen quality in one minutes.
- Large capacity of tray holds a variety of sizes and shapes of containers such as tubes and vials.
- Rubber dimple mat of tray firmly but gently holds containers.
- Test containers can be easily loaded and unloaded without turning the unit off.
- Quality certificate: ISO 9001, ISO 13485.

Specifications:

Model	Fixed speed	Movement	Tilt	Switch	Tray dimension (mm) (WxD)	Weight (kg) (N.W./G.W.)	Overall dimension (mm) (WxDxH)
SR-100	24rpm at 60Hz 20rpm at 50Hz	3D gyratory movement on a central point	20°	ON/OFF	230x150	1.8 / 2.1	230x125x175

Power supply: AC 110V 60Hz (220V 50/60Hz available); The flat mat is available for option.